


Experiment Builder Tutorial

To programme a simple experiment (like the one I showed you) please follow these steps:

1. Click “File” then “New”

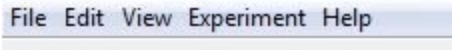

2. Enter a Project Name and select a file location for your experiment then click OK

3. Add any images (or sounds or videos) you want to use to library


- Click “Edit” then “Library Manager”
- Click the + button under the tab for the type of stimuli you want (e.g. under the image tab)
- Select the files you want
- Close the dialogue box

4. Create a display screen for Eyelink set up/ initial instructions:

- Click and drag the **Display Screen** icon from the Action menu to under the start icon



- Click and drag an arrow from the **Start** icon to the **Display Screen** icon

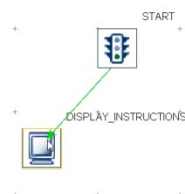


- Double click on the **Display Screen** icon
- Click Insert Multiline Resource (4th button across on the bar at the top)



- Type any instructions you want. Click “Center Align” and “Center Align Vertically” to centre the text in the screen.
- Close the Multiline Resource box
- Right click the Display Screen tab at the top to close, or click the tab that says Experiment

5. Create a keyboard input trigger


- Drag the **Keyboard Trigger** Icon from the trigger menu

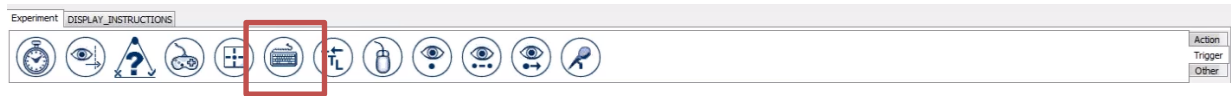

- Click and drag an arrow from the **Display Screen** icon to the **Keyboard Trigger** icon
- To make participants choose only one/specific keyboard key(s):
 1. Click on **Keyboard Trigger** icon
 2. Double click on the box next to where it says “Keys” (should currently say [any])
 3. Select which key you want (e.g. Enter) then close dialogue box
 4. It should now say “[Enter]” next to “Keys”

6. Create Eyelink setup break

- Drag **Camera Setup** icon from the Action menu to under the **Keyboard Trigger** icon

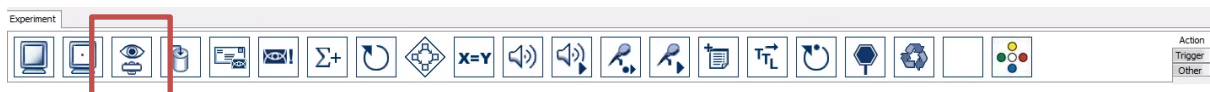

- Click and drag an arrow from the **Keyboard Trigger** icon to the **Camera Setup** icon

7. Create a display screen for experimental instructions

- Repeat step 3, but drag arrow from the **Camera Setup** icon (step 5)

8. Create keyboard input trigger


- Repeat step 4, but drag arrow from second **Display Screen** icon (step 6)

9. Create recorded sequence

- Drag **Sequence** icon from the Action menu, and drag arrow from **Display Screen**



- Click the purple box that appears
- ! Check “Record” in the left hand menu to make sure the Eyelink records the participant’s eye movements and pupil diameter in response to anything in this box. **Box must be ticked to record any eye movements.**

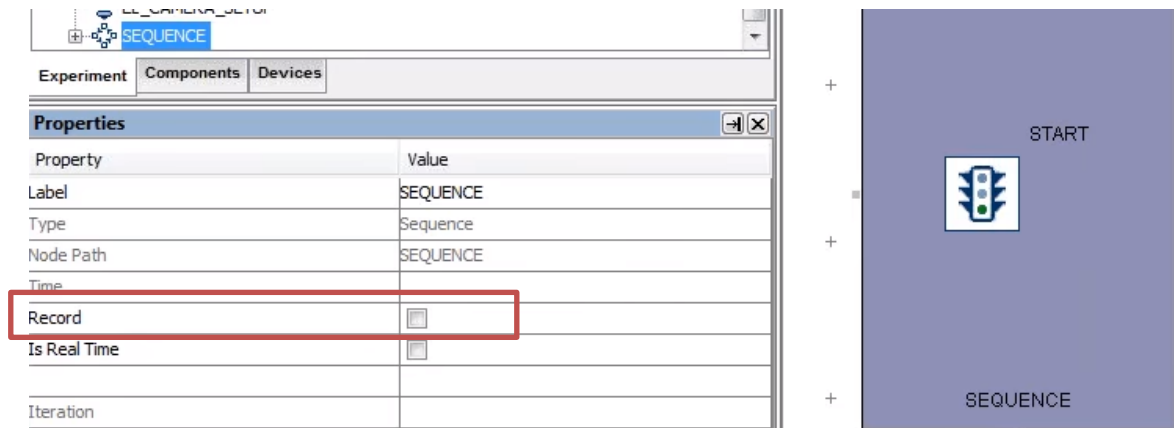

- Double click the purple box to open. This will look like a new experiment, with a new start icon to connect icons to.

In this new sequence:

1. Create Display Screen:

- Drag **Display Screen** from the Action menu
- Connect to the **Start** icon
- Double click the **Display Screen** icon
- Click “Insert Image Resource” (1st button on the bar at the top)



- Select the image you want to insert from the library
 - To centre the image, click on it then click “Horizontal Center Alignment” and then “Vertical Center Alignment” on the bar at the top
 - Right click the Display Screen tab to close
- #### 2. Edit how long the image should stay on screen
- Drag **Timer Trigger** icon from the Trigger menu to under **Display Screen**, and drag arrow from **Display Screen** to connect.
 - Click on Timer icon and change “Duration” to the number of milliseconds you would prefer
- #### 3. Right click sequence tab at the top and close

10. Create screen to thank participants and how to close the screen / finish the experiment

- Repeat step 3, but drag arrow from the purple sequence box

11. Create keyboard trigger to close the experiment

- Repeat step 4, but drag arrow from the step 9 Display Screen

12. Click “Arrange Layout” to make it look nice and neat




Basic Randomisation *(basic but still a little bit complicated!)*

1. Edit the Data Source for the Sequence


- Click the purple **Sequence** box then **Data Source**
- Click “Add Column” then choose a “Column Name” (can be anything, “Stimuli”, “Images” etc.) then choose “String” for “Column Type”
- Add the same number of rows as you have images (3 images=3 rows)
- Type the names of the images into the empty boxes e.g. “Image1.jpg” “Image2.jpg”, exactly as they are in the library
- Click “Randomization Settings”
- Check “Enable Runtime Randomization” and select your column name (e.g. “Images”) for Run Length Control. Then click OK.

2. Create an image that randomly changes

- Double click the **Sequence** box to open it, then double click the **Display Screen** with the image you entered earlier (Step 1 within the sequence box in Step 9)
- Click on the image in the display screen
- In the menu on the left hand side, click “Source File Name”

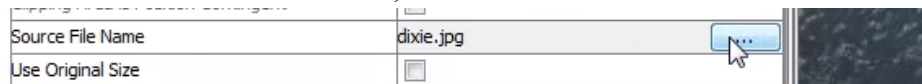

- Then scroll up to **DataSource**

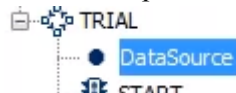

- After selecting this, click the name of the column you created, and then click OK
- Then “Source File Name” should then change to something like (if the columns name was “Image”) @parent.parent.SEQUENCE_DataSource.Image@